


# LASER IMAGER DRYPRO ΣΙΙ

#### Specification of DRYPRO SIGMA 2

| Exposure source                      | Semiconductor laser  |
|--------------------------------------|--|
| Processing method                    | Thermal development system   |
| Film                                 | Dry-imaging recording film SD-S<br>125 sheets per cartridge  |
| Processing capacity                  | $14" \times 17"$ sheets: approximately 70 sheets/hour $8" \times 10"$ sheets: approximately 110 sheets/hour  |
| Pixel size                           | 50 μm (508 dpi)  |
| Network connection                   | Ethernet 10 BASE-T/100BASE-T/1000BASE-T  |
| Film supply                          | Daylight setting   |
| Number of supply trays               | 2ch (Any size 14" × 17", 11" × 14", 10" × 12", 8" × 10") are available   |
| Film size                            | Selectable: 14" $\times$ 17" ( 35 $\times$ 43 cm), 11" $\times$ 14" ( 28 $\times$ 35 cm), 10" $\times$ 12" ( 25 $\times$ 30 cm), or 8" $\times$ 10" ( 20 $\times$ 25 cm) |
| Density adjustment                   | Available  |
| Output grades                        | 14 bit (16384 grades)  |
| Trimmed frame                        | Available  |
| Border                               | Black/Clear  |
| Network connectivity                 | DICOM Print Management SCP   |
| Power                                | 90-130 V (50/60 Hz), 9 A<br>180-264 V (50/60 Hz), 4.5 A  |
| Dimensions (W $\times$ D $\times$ H) | 626 × 649 × 666 mm ( 24.6" × 25.5" × 26.2")  |
| Weight                               | 79 kg  |
| Operating environment                | Temperature: 15-33°C<br>Humidity: 20-80% RH, non-condensing  |

#### Dimensions

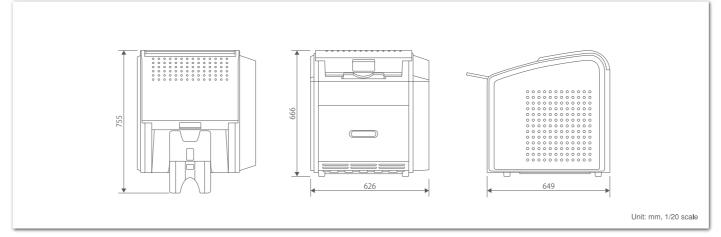

Distributed by:

Brand Name: DRYPRO ∑ II

Type: LASER IMAGER

Model: DRYPRO SIGMA 2

★To facilitate performance improvement, the above specifications are subject to change without notice.



KONICA MINOLTA, INC.

KONICA MINOLTA 1 Sakura-machi, Hino-shi, Tokyo, 191-8511, Japan



LASER IMAGER



160



Konica Minolta has satisfied the strong demand of the medical imaging industry, as the new "DRYPRO ΣΙΙ" dry imager provides optimal diagnostic images.

DRYPRO ΣII creates images of unrivalled clarity and sharpness, utilizing the latest precision optics from Konica Minolta to produce a 50- $\mu$ m pixel pitch;

the highest available resolution.

High Image Quality & Compact Laser Imager

### High-quality images

A semiconductor laser is used to produce a 50- $\mu$ m pixel pitch, resulting in high image resolution. Powerful image-processing algorithms are utilized to simultaneously optimize both the image smoothness and the text sharpness. Diagnostic clarity is preserved and patient data are always legible, regardless of the print size.



## Stable finished quality

The DRYPRO ΣΙΙ density control function maintains the output density via automatic measurement. The system also automatically calibrates whenever a tray of film with a new lot-number is loaded.



## User-friendly design

DRYPRO ΣΙΙ is easy to operate and offers intuitive usability. Film exchange requires simple cartridge insertion. The innovative cartridges allow easy film-size adjustment and support various modalities.





imager. It features a compact footprint as small as  $65 \times 63$  cm<sup>2</sup>, and is designed to integrate into various hospital and clinic.



Four DRYPRO ΣΙΙ film sizes are available:  $14 \times 17$ ,  $11 \times 14$ ,  $10 \times 12$ , and  $8 \times 10$ . Two film sizes can be used simultaneously, for example,  $14 \times 17$  and  $11 \times 14$ , or  $14 \times 17$  and  $8 \times 10$ .



Storing and handling unused film Following confirmation that the film is packed, unused film (like ordinary film) should be stored in a cool, dark place (25°C or below), where it will not be affected by radiation.





